

CE

putting smile by innovative creations

Bionova Corporation, an exclusive global distributor of medical grade Medigel[®]: Organic C₃I Gel Technology based rehabilitation products, Foot Care, Health Care, Gel Cushioning & Operation Theater Product range. The products are exclusively manufactured for Bionova Corporation by Gel Craft Healthcare (P) Ltd.

Technology that raises your lifestyle standards with extreme health benefits.....

C₃I (Cushioning, Comfort and Contour with extreme Inertness) gel with its significant impact resistance, shock absorption (1600 psi); extreme durability; non-toxic; hypoallergenic and doesn't support microbial growth; 100% pressure dispersion ability; high degree of elastomerism (up to 1000% i.e., 10 times relaxed length at failure); pressure dispersion ability; effective pressure relief; fatigue reduction; complete bed sore prevention; bulk specific gravity (0.3 g/cc to 0.9 g/cc) and bio-degradability makes it a versatile technology. At present, the company has commercially available 07 variants of gel technology after~decade old R&D activity with~17500 trials of 285 stabilizers along with 67 quality organic oils. The product range is now available as part of 22 families each one specifies for particular benefits for lifestyle health care.

The technology is developed to provide comfort to all ages during various stages of life and preventive health care benefits on long term usage.

All the devices are developed at WHO, GMP certified manufacturing facility (ISO 9001:2016; ISO 13485:2016). The product ranges are CE complied and accepted by DGCI (India) for commercial applications.



Every step with mediqel® adds calmness to your soul

Medigel[®] **Foot Care Devices** (developed from C_3I gel technology) are extremely soft, comfortable & regained its shape after each step we take. While walking, it conforms to foot shape/bonny prominences & prevents from falling due to non-slip grip.

Its elastomeric property increases the blood flow which reduces blisters and provides relief from various foot pains / ailments. These odorless devices absorb shock, vibrations and prevent frictional impact while distributing body weight evenly. The devices are designed to promote natural healing process & maintain ambient temperature range.

The robust designs are extremely durable that does't tear off or bottom out, an ideal choice for safeguarding feet. It can be easily inserted into any shoe, sandals & foot wears as invisible gel pad, being transparent. These can be easily cleaned with mild shampoo & use talcum powder afterwards.

Medigel® Health Care Devices are designed for pre-hospital trauma care (sealants and bandages); post-surgical care (exercisers, breast prosthesis, scar reduction sheets & knee, elbow & patella protection device); sitting cushions & pads; bed sore sheets & mattresses overlays, infant/ baby care and; interior care (kitchen & garden mats). New generation organic perfumes and body butters (with no synthetic preservative) have recently added by Products.

Percutors personal care. Medipure™ based ultrasound, ECG gels along with microwavable hot and cold packs are launched as new aqua-gel line of products.



Marval Heel Pad

Heel pain reliever gel



pai	ir/	pa	ck

Item	Size	Foot Size
HP-SD	Standard	Universal
HPW-SD	Standard	Universal

- Targets heel pain attributed towards multiple ailments like Plantar Fasciitis or otherwise due to tiredness
- Wedge design elevates the heel for pain relief associated with Achilles Tendonitis
- Heel Pad (wedge) is designed for ladies & children exclusively

Indications

- Sores & Aching Heel Pains
- Plantar Fasciitis & foot pains

Marval Heel Pad (Pando)

Extra cushioning to bony prominences

Item	Size	Foot Size
HN-SD	Standard	Universal









Elastomeric


















- Chronic heel pains (pointing in middle of heel)
- Plantar Fasciitis & Injuries

prominences

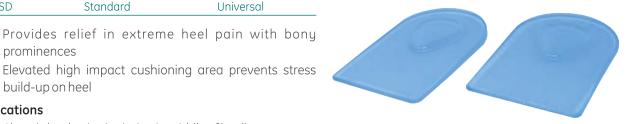

pair/pack

Marval Heel Pad (Donut)

Air cushioning pain reliever gel



- 1	าต	ir/	'n	a	ck

Item	Size	Foot Size
HO-SD	Standard	Universal

- Donut cut out area provides gel aided air cushioning to heel spurs
- Prevents shock absorption by ~ 50% and loss of body alignment

- Pains with Heel spurs (bone outgrowth beneath heel)
- Sensitive Heels





Marval Heel Cup

Cushioning pain reliever gel



Item	Size	Foot Size
HW-M	Medium	Foot Size : 7-8
HW-L	Large	Foot Size : 9-10

- Provides extreme cushioning to heel for relieving heel pain, fat pad atrophy, sores, Plantar Fasciitis and Achilles tendinitis
- Medical grade elastomeric cup contoured to cradle heel cracks & ankle together, thereby prevents side-to-side motion stabilizing foot as it strikes on ground

Indications

- Crack heels and sports injuries
- Achilles tendinitis & arthritic pains

Marval Ladies Heel Cushion

High heel pain reliever

Item	Size	Foot Size
LH-SD	Standard	Universal

- Cushioning to ladies heel while wearing high heels,
- Invisible gel cushion with tapered edges fits easily and alleviates annoying heel pains

Indications

- Fancy shoes / sandals related heel pains
- Plantar Fasciitis and arthritic pains













Impact Stabilizer















pair/pack

Maximal Metatarsal Cushion/Pad

Metatarsal pain reliever





pair/pack

MC-M Medium Foot Size: 7-8 Foot Size: 9-10 Large Supports metatarsal -fore foot area while body

Foot Size

Foot Size: 5-6

- weight is transferred from heel-via-arch to fore foot (gait cucle)
- Useful in management of post fore foot surgery; metastarsalgia and; progressive pain towards ankle

Item

MC-S

Post fore foot surgery & metastarsalgia

Size

Small

Progressive pain towards ankle



Magnatic Arch Support



Item	Size	Foot Size
AR-ES	Extra Small	Foot Size : Junior
AR-S	Small	Foot Size : 5-6
AR-M	Medium	Foot Size : 7-8
AR-L	Large	Foot Size : 9-10

- Realigns body posture by stabilizing foot in neutral position, preventing pronation by contouring with arch plate in progressive flat footers
- Control undesired motion aligning heels & ankles, relieves progressive pain in heel, ankle, knee and spine

Indications

- Inversion / Aversion of Progressive Flat Foot
- Plantar Fasciitis, Achilles tendinitis & arthritic pains

Item	Size	Foot Size
HA-S	Small	Foot Size : 5-6
HA-M	Medium	Foot Size : 7-8
HA-L	Large	Foot Size : 9-10
HA-EL	Extra Large	Foot Size :11-12





Pain Mollifer

















Anti-Slip Grip







- Starter arch with front wedging prevents progressive flat foot, calf muscle stiffness, painful shin splints & progressive pain towards ankles
- Provides posture correctness controls undesired foot motion and prevents transferring of pressure on knees, hip and spine

Indications

- Heel pain with progressive flat foot
- Calf muscle tightening and painful shin splints
- Progressive pain towards ankle

pair/pack

Mollify Insole with Arch

Comprehensive foot care



Item	Size	Foot Size
IA-S	Small	Foot Size : 5-6
IA-M	Medium	Foot Size : 7-8
IA-L	Large	Foot Size : 9-10

- Conforms to foot contours (arch) & supports heel, arch & ankle preventing progressive flat foot, knee and spinal pain
- Enhance body performance by elastomeric blood flow with complete body alignment / body weight re-distribution

- Progressive Flat Foot & metatarsalgic pains
- Achilles tendinitis, Plantar Fasciitis, Foot pains
- Knee/spinal or back/shoulder pains (Body alignment pains)







Mollify Universal Insole

Daily foot insert



item	Size	FOOT Size	
UI-SD	Standard	Universal	

- Thin insole lining conforms to foot shape aid cushioning and prevent foot blisters / sores during long walks / jogs
- Useful for providing relief to all types of foot pain & day long cushioning in all types of shoes / footwear except diabetics and body posture related back pains

Indications

- Lifestyle preventive care for foot pains
- Daily shoe inserts for all type of occupations
- Not useful in foot pains related to body alignment, arthritis or diabetes

Mollify Sports Insole

Foot care while running



- Prevents sports related frictional injuries as starter arch provides median height that makes it ideal for sports and running
- Posture correctness, elastromeric blood flow maintains neutral pronation and impact resistant towards injuries

Indication

- Sports related foot injuries / pains
- Strenuous exercise related foot pains













Impact Stabilizer













Mollify Impact Resistant Insole

Body weight equalizer



Item	Size	Foot Size
Impact Resistant Insole		
IRI-S	Small	5-6
IRI-M	Medium	7-8
IRI-L	Large	9-10
IRI-EL	Extra Large	11-12

- Transfers body weight while equalizing pressure in all directions, elastomeric blood flow, reducing tiredness and pain
- Excellent shock absorber with cushioning to frictional impact of walking & running, safeguards feet against injuries
- Impact Resistant Insole always comes with top fabric layer

Indications

- Daily / intermittent foot pains in middle aged adults, long standing occupational professionals
- Obesity/age related foot, ankle/knee pains

pair/pack



Mollify Diabetic Insole

Neuropathic (sensitive) foot care



pair/pack



- Extreme cushioning for sensitive neuropathic feet, transfers body weight from heel to toe, provides comfort and protection
- Hypoallergenic, temperature homogeneity, antimicrobial layering, impact resistance aids protection of diabetic feet

Indication

- Neuropathic sensitive diabetic foot care
- Knee/spinal or back and arthritic pains

Miracle Diabetic Gel Socks Neuropathic foot care













Impact Stabilizer











- Airflow channel breather design; comfortable soft ribs antiperspiring / microbial fabric with supper soft top (cuff comfort) & soft welt for better grip (non-slip wear) -Elastomeric Oil Dispersion Reflexon Seamless Stitch (EODRSS) Technology.
- Natural touch Medigel® based massaging prevents abrasion, shear forces, callus formation, foot ulcerations, neuropathy of diabetic foot, fat pad atrophy & arthritic pained foot.

- Neuropathic sensitive diabetic feet daily care
- Ideal for old and flabby skin with ankle care



pair/pack

Maximal Toe Separator



Item	Size	Foot Size	
TS-SD	Standard	Universal	

- Cushioning relief from the pain & soreness associated with soft corns or toes rubbing / chafing by reducing shear and interdigital pressure
- Absorbs pressure & friction to end pinching & cramping; prevent aggravation of toe infection due to touch however not recommended to worn inside shoes

- Ingrown toenails protection, corn / blister on toes
- Pinched toes & prevention of infection in diabetic feet



Sk. wali

Maximal Toe Spreader

Straighten & align toes



ոո	ir/	'na	ck

- Item
 Size
 Foot Size

 TR-SD
 Standard
 Universal
- Straighten & align toes by gentle constant pressure while spreading toes, relieving pain of overlapping toes / toe drift and painful bunions.
- Absorbs pressure & friction to end pinching & cramping; prevent aggravation of toe infection due to continuous touch however not recommended to wear inside the shoes

Indication

- Inter digital lesion, big toe pain & overlapping / rubbing toes
- Helomas, blisters, post fore foot surgery and diabetic toes

Maximal Toe Spreader (with D-ring)

Item	Size	Foot Size	
TD-SD	Standard	Universal	

- Pain relieing toe aligner and straightener with extreme cushioning between toes able to provide care similar to Toe Spreader
- Unlike Toe Spreader, its D-ring covers the big toe and thus it can be even used with slippers also. D-ring also provides extra stability to big toe too.

Indications

- Interdigital lesion, big toe pain & overlapping / rubbing toes
- Helomas, blisters, post fore foot surgery and diabetic toes

Big toe pain reliever








































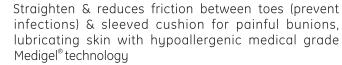

Unlike Toe Spreader, its D-ring covers the big toe and thus it can be even used with slippers also. D-ring also provides extra stability to big toe too.



- Useful for any footwear
- Cushioned sleeve for bunion in shoes









Impact Resistant Gel Socks

Protection against stress fractures



Item	Size	Foot Size
IR-GS-M	M-Medium	Foot Size : 7-8
IR-GS-L	L-Large	Foot Size : 9-10
IR-GS-EL	EL-Extra Large	Foot Size: 11-12

- Designed exclusively for Force Personnel reduces long standing, walking stress while protects feet by fine hand linked toe & cushioned y-heel (Anti-tear resistance) design.
- Gel bed diffuses frictional pressure, impact of uneven land surfaces, prevention of stress (high intense training / exercise) and anti-microbial property aid prolonged wearing during missions.

Indications

- Persons undergoing hard training/exercise
- Ideal for daily use especially for gyms

pair/pack

Miracle Heel Care Socks

Elastomeric protection to high heel wearers



















Easy to Clean



- Elastomeric Oil Dispersion Reflexon (EODR) technology based shock absorbing impact resistant heel protector against high heel shear forces experienced by ladies
- Contoured to provide hind foot grip, heel liner cushioning and nvisible complete heel coverage for high heel wearers

Indications

- Exclusively designed for high heel wearers (sandals or shoes)
- Soften skin reducing sores, cracked heels, blisters, aching heel pains progressive towards ankle, ideal for daily wear

Miracle Sports Gel Socks

Shock absorber against impact injuries



pair	·/pc	ıck



pair/pack

- Item
 Size
 Foot Size

 SS-S
 Small
 Foot Size : 5-6

 SS-M
 Medium
 Foot Size : 7-8

 SS-L
 Large
 Foot Size : 9-10

 SS-EL
 Extra Large
 Foot Size : 11-12
- Designed exclusively for Sports Personnel reduces running stress while protects feet by fine hand linked toe & cushioned y-heel (Antitear resistance) design.
- Reinforced cushioned Y Heel (anti-slip), anti-bacterial / perspiring & reinforced fine hand linked toe for anti-tear resistance (prevent callus formation), available in white color only.

- Sports persons & athletes; exercise induced foot Injuries / joggers
- Non-scheduled game players of varied ages, preventing sports injuries



Flexible Gel Seal

Compression aided care

Coming Soon

Item	Size	Dimensions
FS-SD	Standard	18" x 6" x 5 mm

- Compression gel device with highly absorbing layer on one side while end points-fracture aid bandage & self-sticking velcros able to seal the oozing wound
- Sterile sealant device seals bleeding wounds while stopping blood flow in fields (pre-hospital care)

Indications

- Police, Highways, Army & Paramilitary Patrolling Vehicle
- Ambulances and Evacuation Vehicles
- Emergency Reception (Health Care Units)

pc/pack

Item	Size	Dimensions
WS-SD	Standard	18" x 6" x 5 mm

- Next Generation Medigel® Technology shield with granular coated absorbing layer (on one side) with fracture-aided bandage, supported by velcros to stabilize extreme blood flow
- Radiation sterile shield promotes natural coagulation process by concentrating tissue coagulation factors (removing water from oozing wound)

Indications

- Police, Highways, Army & Paramilitary Patrolling Vehicle
- Ambulances and Evacuation Vehicles
- Emergency Reception (Health Care Units)

Wound Gel Shield





















Posture Erector









Flexible Gel Chest Seal Pre hospital care

 Item
 Size
 Dimensions

 CS-SD
 Standard
 18" x 12" x 5 mm



- Next Generation Medigel® Technology sealant (twice stretchable) with blood absorbing layer stabilizing oozing chest using compression bandage to roll around the chest
- Irradiated field compatible sealant needs to be wrapped around chest like vest & tighten using self-sticking Velcro ends to prepare victim for immediate transfer for definitive care

Indications

- Police, Highways, Army & Paramilitary Patrolling Vehicle
- Ambulances and Evacuation Vehicles
- Emergency Reception (Health Care Units)

pc/pack

Trauma Gel Shield

Stretcher Overlau



- Non-sticky gel sheet for stretcher/field provides cushioning interface for wounded limb while providing pre-hospital trauma care in field/emergency departments
- Natural skin touch, reusable, water proof, extremely soft with pipe corrugation (enhanced surface area) provides ideal overlay for trauma care

Indication

- Highway / Police/ Emergency Patrol Vehicles
- Ambulance Stretchers, ICU/OT Tables



pc/pack

Item	Size	Dimensions
BR-SD	Standard	24" x 36"













Impact Stabilizer





Non Allergenic

- Non-sticky burn roll gel provides cushioning interface to cover the burn site while easily removable without any impact on burnt limb
- Protective non-medicated cover for burnt limb, natural skin touch, water proof and reusable after disinfection using butadiene or similar anti-microbial solution (sterility has to be maintained under medical supervision only)

Indications

- Burn Wards
- Emergency Vehicles/Forces on Special Operations

Flexible Gel Burn Roll Non-abrasive burn care



roll/pack

Flexible Gel Toe Cap




Item	Size	Foot Size
TC-SD	Standard	Universal

- Medical grade cap that isolates the injured/infected toe preventing the aggravation of injury due to frictional impact of the next toe
- Absorbs pressure and friction providing all day relief and protection

- Infected, injured, ingrown toe nail and pinched toes
- Bursa and Lesions on or in between toes

Intense Hand Exerciser

Stress relieving gel ball

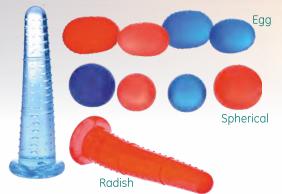

Item	Type/ Size (01 pc/pack)	Design/ Pack
SGB-R-SD	Soft, Standard	Radish
HGB-R-SD	Hard, Standard	Radish
SGB-E-M	Soft/Medium	Egg
SGB-E-L	Soft/Large	Egg
HGB-E-M	Hard/Medium	Egg
HGB-E-L	Hard/Large	Egg
SGB-S-M	Soft/Medium	Spherical
SGB-S-L	Soft/Large	Spherical
HGB-S-M	Hard/ Medium	Spherical
HGB-S-L	Hard/ Large	Spherical

- Enhances blood circulation to extremities increase dexterity and endurance, aids fracture recovery
- 10-15 min daily usage provides positive impact of yog mudras i.e., mritasanjeevani; rudra; prana; chin and vayu
- Aged person and women are advised to use soft types (pink colored) in the beginning; Patients undergoing physiotherapy or fracture recovery should follow radish, egg and then spherical shape gel ball with soft (pink) followed by hard (blue colored) gel balls.

Indications

- Fore-arm fracture / sprained hand recovery during physiotherapy of sportsmen, runners, exercise deprived personnel, diabetic, arthritic, BP patients and students lacking focus
- Uniform personnel or office person especially in middle ages, shaking hands, unusual tiredness, golfer / surgeons / engineers and aged persons with dormant lifestyle.



Bra Size	Dimensions
[30-46; Cup size: A,	(Length *Width
AA, B,C, D,DD]	* Height - in cms)
30A	11.5 x 9 x 3.5
32A 34AA	12.5 x 10 x 4.0
32B, 34A, 36AA	14.0 x 10.5 x 4.5
32C, 34B, 36A, 38AA	14.5 x 11.0 x 5.0
32D,34C, 36B, 38A, 40AA	15.0 x 11.0 x 5.5
32DD, 34D, 36C, 38B, 40A	15.5 x 12.0 x 6.0
34DD, 36D, 38C, 40B, 42A	16.0 x 12.5 x 6.0
36DD, 38D, 40C, 42B, 44A	16.5 x 12.5 x 6.0
38DD, 40D, 42C, 44B	17.0 x 12.5 x 6.5
40DD, 42D, 44C, 46A	17.5 x 12.5 x 7.0
42DD, 44D, 46B	17.5 x 13.0 x 7.0
44DD, 46C	18.0 x 13.5 x 7.5
	[30-46; Cup size: A, AA, B,C, D,DD] 30A 32A 34AA 32B, 34A, 36AA 32C, 34B, 36A, 38AA 32D,34C, 36B, 38A, 40AA 32DD, 34D, 36C, 38B, 40A 34DD, 36D, 38C, 40B, 42A 36DD, 38D, 40C, 42B, 44A 38DD, 40D, 42C, 44B 40DD, 42D, 44C, 46A 42DD, 44D, 46B

Natural Touch Gel Breast Prosthesis





Post-mastectomy care









pc/pack

- Natural skin touch/ skin colored provides elastomeric oil massaging aid in protecting surgically traumatized
- Water resistant (not to remove in swimming pools), no special cradle required to retain shape as impact resistant.
- Hygienic, washable, non-leaky (as other jelly breasts) with long life utility.

- Aid in providing artificial breast during post-mastectomy (surgical removal of complete breasts)
- Should not be used if artificial breast implant is in place as it is only external-use device







Magic Scar Reduction Sheet



Item	Size [Length *Width * in inches]	Available Packs 1, 2, 5, 12 pc/pack
SRS-12	1" x 2"	2X, 5X, 12X
SRS-18	1" x 3"	1X, 2X, 5X
SRS-25-1	1" x 4"	1X, 2X,5X
SRS-32	1" x 5"	1X, 2X,5X
SRS-38-1	1" x 6"	1X, 2X,5X
SRS-20	1.5" x 2"	1X, 2X,5X
SRS-30	1.5" x 3"	1X, 2X,5X
SRS-40	1.5" x 4"	1X, 2X,5X
SRS-50-1	1.5" x 5"	1X, 2X
SRS-60	1.5" x 6"	1X, 2X
SRS-25-2	2" x 2"	1X, 2X
SRS-38-2	2" x 3"	1X, 2X
SRS-50-2	2" x 4"	1X, 2X
SRS-62	2" x 5"	1X, 2X
SRS-75-1	2" x 6"	1X, 2X
SRS-75-2	4" x 3"	1X, 2X
SRS-100	4" x 4"	1X, 2X
SRS-125	4" x 5"	1X
SRS-150	4" x 6"	1X
SRS-112	6" x 3"	1X
SRS-188	6" x 5"	1X
SRS-225	6" x 6"	1X

1" ~2.5 cm; 1.5" ~4cm; 2"~5cm; 4"~10 cm & 6"~10 cm. (1X- 1pc/ pack, 2X- 2 pc/pack, 5X- 5pc/pack, 12 pc/pack)

Elastomeric







Indications



Non Allergenic

Tummy tuck (abdominoplasty) & Buttock augmentation

providing inert and hygienic micro-environment

C-section & Stenotomy scars; brachioplasty; breast surgery

Post surgery hypertrophic and keloid scars on any anatomical body part

Organic oil moistening reduces irritation/ itching, promotes natural healing & faster recovery by

Elastomeric compression aids in reducing hypertrophic and keloid scars by flattening, softening,

Easily removable, soothing natural skin touch and leaves no greasy marks (non-fluidic). Use it as per

smoothing & bending of scarred tissue in 8-12 weeks (depending upon age & type of scar)

instructions given in the instruction booklet in the pack. It should not be used on open wounds

CEM Medigel Care

Gel Knee Ring

Knee Protection Device Insert





pc/pack

Item	Size	Pack Size
OEM-PS-CS	Customizable	OEM Supply

- C shape cut out designed to support patella from hind side to provide shock resistant.
- Gel protection device against short time high impact shock.
- Various dynamic variants are available based on customer's requirement for OEM supply only.

Indication

Aid is patella support / knee protection device

Tennis Elbow Pad

Type/ Size Pack Size Item Oval/Standard **OEM Supply KRO-SD** KRR-SD Round/Standard **OEM Supply**

- Medigel® Knee Ring is designed to get fit in Knee Support devices.
- Support to knee cap-open patella & protects knee by restricting unwanted movements.
- Aid in body alignment & cushioning protective interface is provided by the Medigel®.

Indication

Aid in relieving from Back Pains / Knee Pains

Pattela Support

Pattela Protection Device Insert



Pack Size

OEM Supply



pc/pack

Impact Stabilizer

Anti-Slip Grip

Posture Erector









Medigel® device insert for fitting into Tennis Elbow

- Support Device.
- Aid in relieving the pain by putting right strain while cushioning elbow area simultaneously

Indication

Item

TEP-SD

Relieving from Tennis Elbow Pains

Size

Standard



Sorry to hear, someone broke his rear, it's pain one need to bear, it requires special care while the one is sitting on the chair



Chillax Coccyx Seat Cushion

Tail injury care

Item	Type/Configuration/Top Gel*	Dimensions
	layer~18 mm);	L x B x Tf thickness in front) x Tb (Thickness at back) {in inches}]
CSC-F-C-SD	Flat/Corrugated/Standard	16" x 16" x 2.5" x 2.5
CSC-W-C-SD	Wedge/Corrugated/Standard	14" x 16" x 3.75" x 2

Features

- Seat cushion with coccyx cut out to provide air cushioning to tail bone, top corrugated sheet distribute pressure equally at buttocks, manages body heat and conforms body shape.
- Sitting mass, stability to cushion is provided by closed cell high density PU foam beneath corrugated gel sheet. encapsulated in incontinent-proof and breathable extenson fabric (wicks away moisture & heat)
- Wedge type (with slope of 1.75") allows weight transfer to feet faster, aid spinal alignment, though if your pelvic structure or prolonged sitting seat is flat, flat type is advisable.

Indications

- Designed exclusively for coccyx pain relief
- Pipe type configuration (corrugation) allows continuous flow of air and thereby anti-sweating



Chillax Sitting Seat Cushion











Item	Type/Configuration/Top Gel	Dimensions
	layer~18 mm	L X B X Tf thickness in front) X Tb (Thickness at back) {in inches}]
SSC-F-BX-SD	Flat/Box Type/Standard	16" x 16" x 2.75" x 2.75
SSC-F-B-SD	Flat/Ball Type/Standard	16" x 16" x .75" x 2.75"
SSC-W-C-SD	Wedge/Corrugated/Standard	16" x 14" x 2.25" x 1.25"



















Features

- Top corrugated gel layer (20 mm) distribute pressure equally at buttocks, manages body heat, enhanced blood circulation and conforms body shape.
- Sitting mass, stability to cushion is provided by closed cell high density PU foam beneath gel sheet, encapsulated in incontinent-proof and breathable extenson fabric (wicks away moisture & heat)
- Wedge type (with slope of 1") allows weight transfer to feet faster, aid spinal alignment, though if your pelvic structure or prolonged sitting seat is flat, flat type is advisable.

- Reduces prolonged sitting fatigue
- For all types of chairs (Official and Home use)
- Dining Chairs / Car Seats
- Wheel Chair





Chillax Seat Cushion for Cars

Anti-vibrational/fatigue prevention care



pc/pack

Item	Type/Configuration/Top Gel	Dimensions
	layer~18 mm);	L X B X Tf thickness in front) X Tb (Thickness at back) {in inches}]
SSC-W-C-SD	Wedge/Corrugated/Standard	14" x 16" x 2.25" x 1.25"

Features

- Anti vibrational top corrugated gel layer distribute pressure equally, enhanced blood circulation and absorbs 95% torque pressures
- Closed cell high density PU foam base, item is encapsulated in incontinent-proof and breathable extenson fabric (wicks away moisture & heat)
- Wedge type (with slope of 1") allows weight transfer to feet faster, aid spinal alignment,

Indications

Designed exclusively for driving seat of automobiles



Chillax Wheel Chair Cushion

Blood Circulation Enhancer





Item	Type/Configuration/Top Gel	Dimensions
	layer~18 mm	L x B x Tf thickness in front) x Tb (Thickness at back) {in inches}-2.75" x 2.75"]
WCC-F-C-M	Flat/Corrugated/Medium	14" x 14"
WCC-F-BX-SD	Flat/Box Type/Standard	16" x 16"
WCC-F-C-L	Flat/Corrugated/Large	18" x 18"
WCC-F-B-M	Flat/Ball Type/Medium	14" x 14"
WCC-F-B-SD	Flat/Ball Type/Standard	16" x 16"
WCC-F-B-L	Flat/Ball Type/Large	18" x 18"

Features

- Top corrugated gel layer distribute pressure equally at buttocks, manages body heat, enhanced blood circulation and conforms body shape.
- Designed in ball type (bouncing elastomerism); box type (anti-tiring) & corrugated (anti-sweating) with sitting sore prevention medigel technology

Indications

Exclusively designed for Wheel Chairs only

















Chillax Wheel Chair Pads

Blood Circulation Enhancer



Item	Configuration/Gel Layer	Dimensions
	Thickness~18 mm);	Length x Width {in inches}]
WCP-BX-SD	Box Type/Standard	16" x 16"
WCP-B-SD	Ball Type/Standard	16" x 16"
WCP-B-L	Ball Type/Large	18" x 18"
WCP-C-L	Corrugated/Large	18" x 18"

Features

- Gel layer distribute pressure equally at buttocks, manages body heat, enhanced blood circulation and conforms body shape.
- Prevents pressure on pelvic structure during prolonged Wheel Chair sitting

Indications

- Designed exclusively to lay over the wheel chair cushions, available in multiple configurations
- Ball type (bouncing elastomerism); box type (anti-tiring) and; corrugated (anti-sweating) sitting sore prevention medigel technology

Chillax Piles Seat Cushion

Elastomeric donut-aided air cushioning











Impact Stabilizer



Equal Pressure Decribution





Clean Clean

A Non Allergenic

ItemType/Configuration/Top Gel
layer~18 mm);Dimensions
L x B x Tf thickness
in front) x Tb (Thickness at
back) {in inches}-2.75" x 2.75"]PSC-F-C-SDFlat/Corrugated/Standard16" x 16PSC-F-BX-SDFlat/Box Type/Standard16" x 16

Features

- Flat seat cushion with donut cut out to provide air cushioning to middle of buttock area (piles), top corrugated / box type sheet distribute pressure equally at buttocks, manages body heat and conforms body shape.
- Sitting mass, stability to cushion is provided by closed cell high density PU foam beneath corrugated gel sheet, encapsulated in incontinentproof and breathable extenson fabric (wicks away moisture & heat)

Indications

Exclusively designed for piles patients



Extense Gel Overlays (Plain gel sheets)

Cushioning that comforts you



Item	Size/Soft	Dimensions L X B X (in inches)
X thickness (in	mm)	
FGS-S-3-SD	3 mm/Soft	24" x 24 " x 3 mm
FGS-S-6-SD	6 mm/Soft	24" x 24 " x 6 mm
FGS-S-9-SD	9 mm/Soft	24" x 24 " x 9 mm
FGS-S-12-SD	12mm/Soft	24" x 24 " x 12 mm

Features

- Sheet Thickness increases cushioning as 3, 6, 9 and 12 mm can reduce tossing and turning pressure; antivibrational; bed sore prevention and; fatigue /impact resistant respectively and cumulatively
- Pressure point reliever with one sided non-woven sheet

Indications

- Protective cushioning & gellining (used in artificial shoes)
- Multiple uses in wide range of products and devices

Rem Mattresses Overlay

3 pc set /pack









Type/Configuration/Top Layer Gel Thickness (mm)~18mm (3 pieces set of same size)	Dimensions [Length \times Width \times T $_{\rm b}$ (back thickness) \times T $_{\rm f}$ (front thickness) {in inches}~2.7"]
Flat/Corrugated/Standard	36" x 24"
	(3 pieces set of same size)



- 03 piece set (encapsulated in zipped lycra) with anti-sweating corrugation (18 mm) and extremely hygienic/portable
- Extremely durable, simple dip washable and useful on any mattresses

Indications

- Extremely useful in bed sore patients
- Daily use as overlay sheet for home mattresses / field picnics



Best way to sleep



Anti-Slip Grip



Easy to Clean





Rem Bed Sore Sheet

Sleeping health care

Item	Configuration/Size/Gel Thickness (mm) (3 pieces set of same size)	Dimensions [Length x Width {in inches}]
BSS-P-SD-12	Plain/Standard/12mm	36" x 24"
BSS-C-SD-18	Corrugated/Standard/18mm	36" x 24"

Features

- Anti-sweating (corrugated) or extreme pressure distribution (plain) variants of bed sore sheets provide ideal choice to patients suffering or prone to suffer from bed sores
- Extremely softness grade of gel provides enhanced comfort and care to patients

Indications

- Extremely useful in ICU, bed ridden patinets getting hospital/homecare
- Useful technology for comatose patients too



pc/pack





Pressure Relieve Gel Overlays (Customized)

Customize your desires







Impact Stabilizer



Petar



Clean





Features

- Provides extreme cushioning to any traumatized part while taking first aid at emergency section/ field attributed towards complete elastomeric pressure distribution.
- Available in multiple variants for multiple uses especially bed sore prevention

Indications

• Individual and Industrial Multiple uses

Variants-

Softness: Shore A 00 to 05

Hardness: Shore A10-15 (footwears)

Color: Natural-blue, pink or green

Configuration: Plain, ball or corrugated

Style: Non woven of fabric cover, one or both sides or one

side fabric, other side non-woven

Size: Minimum size: 2" X 2" -to- Maximum

size: 24" X 24"

Thickness: 3mm to 20mm (Average)

Children are the hands by which we take hold of heaven Henry Ward Beecher

Cozy Baby Pillow

Protection against stress fractures

Item	Size	Foot Size
BP-SD	Standard	01 pc/pack

- Alignment of soft skull and extreme elastomeric cushioning aids baby's head
- Extremely soft natural skin touch with oil massaging impact

Indications

- Aids in protection of soft skull and undesirable movementsin daily usage
- Conforms with soft neck with complete pressure redistribution



pc/pack

Cozy Baby Gel Underlay

All night father's hug care



Organic























ه. ا



- Developed from bed sore prevention technology with temperature management and unlike rubber baby beds, no slipping of bed sheet
- Water proof technology with friction less urine flowing out of overlay prevents wetting of baby for longer periods during sleep
- Elastomeric pink color baby bed reduces tossing / turning pressure, doesn't bottoms out and one side covered with non-woven layer, simple dip washable


- Aids in protection of baby's body and induces comfortable sleep due to complete pressure redistribution
- Daily usage aids oil massing impact moistened soft skin of baby

Cozy Neonatal Gel Underlay

Natural cushioning care

Item	Dimensions [L" x W"] Thickness (in mm)]	Available Pack
NGU-24	24" x 18" x 24 mm	01 pc/pack

- Developed exclusively for new born babies (infants / neonatal) with ball type configuration provide air cushioning to soft body parts
- Water proof and gaps between balls take urine that flows out prevents wetting of baby for longer periods. Bouncing impact of ball type provides soft massaging impact.
- Elastomeric pink color neonatal bed disperses pressure, doesn't bottom, one side covered with non woven fabric and simply dip washable

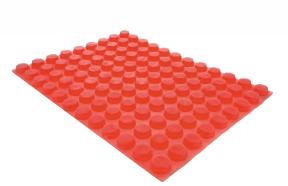


pc/pack

Indications



Aids in sleep of neonatal (infant) by protecting body via dispersing pressure in all directions inducing comfortable sleep



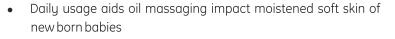





















Item	Size	Pack Size
PR-SD	Standard	10 ml

- Foot Roll on perfumes are developed for Medigel® foot/ health care devices being completely organic with no synthetic additives or alcohol
- Use of any alcohol based perfume is prohibited with Medigel® devices, thus just apply this organic perfume on your foot / Medigel® device

- Use with Medigel® Foot care / health devices
- Personal Foot Care







pc/pack



- Apply fragrance on any area you want to smell good. The inside of your wrist and behind the year works well.
- Used along with gel socks, insoles or any other devices.
- Doesn't contain water, mineral oil, preservatives, artificial dyes, alcohol, fillers, parabens, phthalates and formaldehyde

Item	Size	Pack Size
QS-SD	Standard	6 ml (0.20 fl.oz)
QS-P	Personal	1 ml (0.03 fl.oz)







Anti-Slip Grip

Easy 70 Olean Easy to Clean

Ą Non Allergenic

pc/pack



Once a memory association is created from a scent and is then connected with a person, place and eventit's hard to associate it with something else. Dr. Palemera Dalton, Olfactory Scientist, Philadelphia



Serene & astonishing, just the way you are



Item	Size	Pack Size
CBB-M	Medium	50 gm
CBB-L	Large	170 gm













Impact Stabilizer













- Developed with blend of Kokum and multiple butters & essential oils with no synthetic preservatives, dye, alcohol, fillers, parabens etc.
- Hydrates & soften your skin with day lasting smoothness for entire body.

Indications

- Daily use for softening your skin
- Personal, Spas, Cosmetic clinics, Swimming pools and Yoga centres


Item Size **Pack Size** TO-SD Standard 6 ml (0.20 fl.oz)

- The essence of its first note is beautiful garden of rose with refreshing mint and lavender flowers in the background.
- This formulation is developed to enhance your pheromones to blend with it and give you a unique fragrance.

Indications

- Apply fragrance on any area you want to smell good. The inside of your wrist and behind the year
- Used along with gel socks, insoles or any other devices.
- Doesn't contain water, mineral oil, preservatives, artificial dyes, alcohol, fillers, parabens, phthalates and formal dehyde

Cokum Body Butter

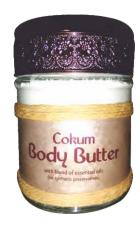

pc/pack

Cokum Foot Cream **Healing Care**




- **Pack Size** Item Size CFC-SD Standard 20 gm
- Organic blend of Kokum, multiple butters & healing essential oils g with no synthetic ingredient / preservative
- Heals your foot, removes unpleasant smell, soften your skin and pain reliever / prevents infections

- Daily use for keeping your foot healthy / relieves pain
- Athletic foot, diabetic, arthritic, old / flabby skin foot
- Foul Smelling foot care



Once a memory association is created from a scent and is then connected with a person, place and eventit's hard to associate it with something else. Dr. Palemera Dalton, Olfactory Scientist, Philadelphia



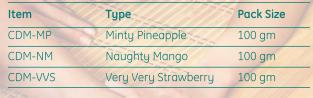

- Restorative non-comodogenic sebum balancer hydrating skin protector
- Non clogging pores with extreme skin care by enriched vitamin E emmollient

Indications

- Daily usage keep your skin soft and supple
- Minty fruity variants tantalize your sensations
- Personal, Spas & problems of dry skin



pc/pack

Item	Size	Pack Size	
SBB-M	Medium	50 gm	
SBB-L	Large	170 gm	




- Gentle healing with non saponifiable components soften skin as it melts at skin temperature
- Comforting & cooling for extremely soft & sensitive skin

Indications

- Relieves inflammation with continued topical usage
- Personal, Spas, Cosmetic Clinics and Yoga Centres


Shea Body Butter Camforting skin care

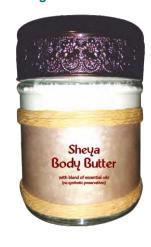

pc/pack

Shea Deep Moisturizer

Cooling skin care

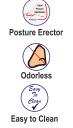

Non Allergenic

Impact Stabilizer

Anti-Slip Grip



Item	Туре	Pack Size
SDM-MP	Minty Pineapple	100 gm
SDM-NM	Naughty Mango	100 gm
SDM-VVS	Very Very Strawberry	100 gm

- Soothes sunburn, heals dry skin & anti-inflammatory
- Reduces eczema, peeling & stimulate collegen production (youth scaffolding protein)

- Cooling massage keep your skin healthy and vibrant
- Minty fruity variants tantalize your sensations
- Personal, Spas & skin healing care



"Good Health is the only blessing we all want, while time is only currency we can't earn back"

medi pure™ Ultrasound gel





pc/pack

Item	Wt	Pack Size
MP-ECG-S	260 g (9.17oz)	25 Pcs /Pack
MP-ECG-L	183.40 oz (5.2 Kg)	6 Pcs / Pack

- High Transparent Gel with optimum viscosity.
- Non corrosive formulation prevents the electrodes from corrosion and increases their life span.
- Water based Gel.
- Ideal for electrodes as it can be washed off easily.

Indications

- Health / Diagnostic Clinics
- Hospital Facilities

Item	Wt	Pack Size
MP-UG-S	260 g (9.17oz)	25 Pcs /Pack
MP-UG-L	183.40 oz (5.2 Kg)	6 Pcs /Pack

- High quality acoustic couplant developed Specially for Sonography
- Non-greasy, Non-Irritant, Non-toxic, Water Soluble
- Easily removable with tissue paper.

Indications

- Physiotherapy / Diagnostic Clinics
- Hospital Facility

medi pure™ Advanced ECG gel

Pack Size

1 Pc/Pack

1 Pc/Pack

1 Pc/Pack

1 Pc/Pack













pc/pack

Elastomeric








Posture Erector







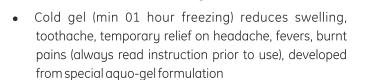

Shoulder Wedge

Wrap Size

Small

Large

Medium

 Hot Gel (30 sec microwave or up to 4 min in boiling water) provides temporary relief of muscle aches and pains, back pains, foot / shoulder pains, sinusitis and menstrual cramps



Item

HCP-S

HCP-M

HCP-L

HCP-SW

- Domestic (home use)
- Physiotherapy & Recovery Clinics





pc/pack

medi pure™ Cold Pack Alleviating Pains



Item	Code Size	Pack Size
CP-SD	Standard	01 Pc /Pack

- Translucent Cold Pack is designed for providing edema related pains; headaches, fever and burning pains
- Highly durable, long lasting and extremely easy to use, aqua-gel formulation achieves freezing temperature quickly

Indications

- Domestic (Home/Personal Use)
- Physiotherapy & Recovery Clinics

pc/pack

Item	Code Size	Pack Size
GM-SD	Standard	01 Pc /Pack

- Provide protection to knees while working in gardens in bending position for long periods
- Safeguards knees while seating on rough surfaces on gardens in bending positions attributing to shock absorption ability (1600 psi)



- Gardeners, Domestic Use, Green Houses, Parks Maintenance
- Automobile Service/ Repair & Masons
- Any occupation which require knee bending position for long time

Garden Kneeler Gel Cushion

Gardening made easy

Pack Size

01 Pc/Pack



pc/pack

Interior Care Kitchen Mat

No more leg pains while cooking



 Kitchen mats are exclusively designed for ladies standing in kitchen for long time

Code Size

Standard

- Encapsulated with gel layer aided cushioning, elastoemeric enhanced blood flow and antitiring impact
- Also useful in events, ticket purchasing queues or security personnel standing for long time at one location

Indications

Item

KM-SD

- Long Standing occupations like catering services; pandal queues, ticket counter fronts, Security Personnel Deployed Locations etc.
- Designed for Modulated Kitchens and Classical Cooking Areas

pc/pack



Easy to Clean

Elastomeric

Impact Stabilizer

Anti-Slip Grip



Organic C₃I Gel Technology vs. Other Technologies

Feature	Organic C₃I Gel	Silicone Gel	Others (EVA & PU Foam)
Definition (touch/leak)	Natural skin touch oil based inert elastomers.	Inert & synthetic gelly based formulations with 4% chances of allergenicity & accidental leaks in long term use	Polymeric ethyl vinyl acetate foam material. Its touch is synthetic away from natural touch however don't leak, with time it loses its air compression capability (PU Foam)
Cushioning & Comfort	Far superior than Silicone Gel	Far superior than rubber based products	Less cushioning & comfort with limited shelf life
Elasticity/ Elastomerism	Extensive elasticity and elastomeric bounce back property	More Elastic than EVA or PU based products, with least elastomerism	EVA has more elasticity than PU Foam however lesser than counterparts, loss of bouncing backwith use
Composition & Development	Organic oils mixed as per specified patented formulation, and finished products are developed by injection and moulding technologies for desired purpose at extremely higher temperature conditions.	Powder form is mixed with other reactants to develop a polymer of desired use and by suitable cast-fit technology, taken into various shapes.	Polymeric Sheets of EVA able to get into various forms through engineering aid is possible while PU Foam has air gaps which undergoes compression and decompression cycle
Combustion	Being elastomer it melts like a liquid with an invisible fume and nearly odorless. The burning gel has significant higher temperature due to phase change and should not be touched.	When silicone is burned in air or oxygen, it forms solid silica (silicon dioxide) as a white powder, char, and various gases. The readily dispersed powder is sometimes called silica fume	Ethyl Vinyl Acetate (EVA) – extremely toxic in nature while Poly Urethane (PU) Foam also produces polluted combustion air thus requires to be used cautiously

Note: Keep off sharp/pointed implements. Do not expose to heat or direct sunlight.

Instructions for Care

- 1. Keep it in cool place and wash it with any detergent based solution. Dip it in the solution for 10-15 minutes, remove the dirt with hand towel, wash it fresh running water.
- 2. Don't use sharp objects like scissors or pins etc.
- 3. Don't use perfumes, sprays especially those containing alcohol.
- 4. Dry the washed items under shade area.



BIONOVA CORPORATION

A UNIT OF GELCRAFT HEALTHCARE PVT. LTD.

WHO GMP Certified Facility

(ISO 9001 : 2016 & ISO 13485 : 2016)

A unit licenced sole Indian Collaborator of Gel Craft (Italy)

MFG Faictility India:

A2/48, Site 4, Sahibabad Industrial Area, Ghaziabad-201010, Uttar Pradesh, India Mobile: +91 9810717192 | Email: gelcrafthealthcare@gmail.com